

A. Midesophageal Views

Acquisition Protocol

Transducer Angle: - 0° - 10°
Level: Mid-esophageal
Maneuver: (from prior image):
NA

1. **ML 5 Chamber View**

Transducer Angle: - 0° - 10°
Level: Mid-esophageal
Maneuver: (from prior image):
Advance & Rotate

2. **ML 6 Chamber View**

Transducer Angle: - 50° - 70°
Level: Mid-esophageal
Maneuver: (from prior image):
NA

3. **ML Mixed Cross-sectional View**

Transducer Angle: - 80° - 100°
Level: Mid-esophageal
Maneuver: (from prior image):
NA

4. **ML 2 Chamber View**

Transducer Angle: - 120° - 140°
Level: Mid-esophageal
Maneuver: (from prior image):
NA

5. **ML Long Axis View**

Transducer Angle: - 120° - 140°
Level: Mid-esophageal
Maneuver: (from prior image):
Withdraw & Anteflex

6. **ML A4 LAX View**

Transducer Angle: - 90° - 110°
Level: Upper-esophageal
Maneuver: (from prior image):
Withdraw

7. **ML Ascending Aorta LAX View**

Transducer Angle: - 0° - 30°
Level: Upper-esophageal
Maneuver: (from prior image):
Clockwise


8. **ML Ascending Aorta SAX View**

Transducer Angle: - 0° - 30°
Level: Upper-esophageal
Maneuver: (from prior image):
CW, Advance

9. **ML Right Pulmonary Vein View**

Transducer Angle: - 25° - 45°
Level: Upper-esophageal
Maneuver: (from prior image):
CCW, Advance, Anteflex

10. **ML A4 SAX View**

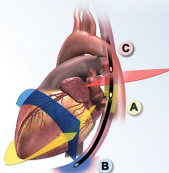

Acquisition Protocol

Transducer Angle: - 50° - 70°
Level: Mid-esophageal
Maneuver: (from prior image):
CW, Advance

11. **ML RV Inflow-Outflow View**

Transducer Angle: - 50° - 70°
Level: Mid-esophageal
Maneuver: (from prior image):
CW

12. **ML Moored Biventricular TV View**

Transducer Angle: - 90° - 110°
Level: Mid-esophageal
Maneuver: (from prior image):
CW

13. **ML Biventricular View**

Transducer Angle: - 90° - 110°
Level: Upper-esophageal
Maneuver: (from prior image):
Withdraw, CCW for the right veins, CCW for the left veins

14. **UL Right and Left Pulmonary Vein View**

Transducer Angle: - 90° - 110°
Level: Mid-esophageal
Maneuver: (from prior image):
Advance

15. **BL Left Atrial Appendage View**

B. Transgastric Views

Acquisition Protocol

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Advance & Anteflex

16. **TG Basal SAX View**

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Advance & Anteflex

17. **TG Mid Papillary SAX View**

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Advance & Anteflex

18. **TG Apical SAX View**

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Anteflex

19. **TG RV Basal View**

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Right-flex

20. **TG RV Inflow-Outflow View**

Transducer Angle: - 0° - 20°
Level: Transgastric
Maneuver: (from prior image):
Left-flex, Advance, Anteflex

21. **Deep T 2 Chamber View**

Transducer Angle: - 90° - 110°
Level: Transgastric
Maneuver: (from prior image):
Neutral flexion, Withdraw

22. **TG 2 Chamber View**

Transducer Angle: - 90° - 110°
Level: Transgastric
Maneuver: (from prior image):
CW

23. **TG RV Inflow View**

Transducer Angle: - 120° - 140°
Level: Transgastric
Maneuver: (from prior image):
CCW

24. **TG LAX View**

C. Aortic Views

Acquisition Protocol

Transducer Angle: - 0° - 10°
Level: Transgastric to mid-esophageal
Maneuver: (from prior image):
Neutral flexion

25. **Descending Aorta SAX View**

Transducer Angle: - 90° - 100°
Level: Transgastric to mid-esophageal
Maneuver: (from prior image):
Neutral flexion

26. **Descending Aorta LAX View**

Transducer Angle: - 0° - 10°
Level: Upper-esophageal
Maneuver: (from prior image):
Withdraw

27. **UL Aortic Arch LAX View**

Transducer Angle: - 75° - 90°
Level: Mid-esophageal
Maneuver: (from prior image):
NA

28. **UL Aortic Arch SAX View**

Mitral Valve Views

Diagram of Mitral Valve Views:

The diagram shows the mitral valve in cross-section with the following views and angles:

- P₁-A₁A₁-P₁** (0°)
- A₁A₁-P₁P₁** (0°)
- P₂-A₂** (45°)
- P₃-A₃** (90°)
- P₄-A₄** (135°)
- P₅-A₅** (180°)
- P₆-A₆** (225°)
- P₇-A₇** (270°)
- P₈-A₈** (315°)

The diagram also shows the mitral valve in cross-section with the following views and angles:

- P₁-A₁A₁-P₁** (0°)
- A₁A₁-P₁P₁** (0°)
- P₂-A₂** (45°)
- P₃-A₃** (90°)
- P₄-A₄** (135°)
- P₅-A₅** (180°)
- P₆-A₆** (225°)
- P₇-A₇** (270°)
- P₈-A₈** (315°)

3D TEE Image Acquisition

Structures Imaged	Acquisition Protocol	Structures Imaged	Acquisition Protocol	Structures Imaged	Acquisition Protocol
Left Ventricle (Acquire using multi-beat mode)	- Obtain a view of the left ventricle from the 0°, 60°, or 120° mid-esophageal positions. - Use the biplane mode to check that the left ventricle is centered in a second view 90° to the original.	Interventricular Septum (Acquire using single-beat or multi-beat mode)	Obtain a view of the interventricular septum from the 0° mid-esophageal with the probe rotated to the interventricular septum.	Mitral Valve (Acquire using multi-beat mode, single-beat mode)	- Obtain a view of the mitral valve from the 0°, 60°, 90°, or 120° mid-esophageal views. - Use the biplane mode to check that the mitral valve annulus is centered in a second view 90° to the original.
Right Ventricle (Acquire using multi-beat mode)	Obtain a view of the right ventricle from the 0° mid-esophageal position with the right ventricle tilted so that it is in the center of the image.	Aortic Valve (Acquire using single-beat or multi-beat mode)	Obtain a view of the aortic valve from either the 60° mid-esophageal, short-axis view or the 120° mid-esophageal, long-axis view.	Tricuspid Valve (Acquire using multi-beat mode, single-beat mode)	Obtain a view of the tricuspid valve from either the 0° to 30° mid-esophageal, 4-chamber view then so that the valve is centered in the imaging plane in the 60° transgastric view with anteflexion.